

IMP-PX-100  
Version F

# OXY-30000

EN



OXYPHARM

829 rue Marcel Paul  
94500 Champigny-sur-Marne  
commercial@oxypharm.net  
T: +33 1 45 18 78 70

[www.oxypharm.net](http://www.oxypharm.net)



I.	<b>INTRODUCTION</b>	2
I.1.	General information	2
I.2.	Classification / Regulation	2
I.3.	Warranty	2
I.3.1.	Clause	2
I.3.2.	Purpose	2
I.3.3.	Transport	2
I.3.4.	Duration	2
I.4.	Technical characteristics	3
II.	<b>PRECAUTIONS FOR USE</b>	3
II.1.	Safety instructions	4
III.	<b>DESCRIPTION</b>	5
III.1.	OXY-30000 description	5
III.2.	OXY-30000 dimensions	5
IV.	<b>INSTRUCTIONS FOR USE</b>	6
IV.1.	Positioning of the device	6
IV.2.	Installation of the bottle	6
IV.3.	Diffusion duration and consumption	7
IV.4.	Functioning in standard mode	8
IV.4.1.	Use without remote control	8
IV.4.2.	Use with remote control	9
IV.4.2.1.	Pairing with Nocospray (« blank » remote control)	9
IV.4.2.2.	Functioning with the remote control	10
IV.4.2.3.	Pairing with another Nocospray (« re-pairing »)	11
IV.5.	Functioning in timer mode	12
IV.5.1.	General informations	12
IV.5.2.	Timer mode features	13
IV.5.3.	Back on standard mode	13
V.	<b>TRACEABILITY SYSTEM</b>	14
V.1.	Software installation	14
V.2.	Transferring treatments historical	15
V.3.	READING AND MODIFYING DATA	16
V.3.1.	General information	16
V.3.2.	General settings modification	17
V.3.3.	Export records (treatment historical)	17
V.3.4.	USB mode exit	17
VI.	<b>TROUBLESHOOTING</b>	18
VII.	<b>HYGIENE AND MAINTENANCE</b>	18
VIII.	<b>EMC COMPLIANCE BOARD</b>	19
IX.	<b>RECYCLING</b>	21

# **I. INTRODUCTION**

## **I.1. General information**

This user manual describes the functioning of the OXY-30000 diffuser and includes essential information to ensure the safe use of the machine.

The documentation applies to Nocospray 2.0 diffusers in their original manufacturer's condition and which have not been altered.

Review this manual before using the OXY-30000 diffuser. It was written to ensure you master all functionalities.

## **I.2. Classification / Regulation**

Electric classification : Class I

Compliant to CEM requirements according to standard 60601-1-2 (2007)

## **I.3. Warranty**

The warranty certificate should be submitted to the dealer or to OXY'PHARM company, Z.A des Grands Godets, 829 rue Marcel Paul, 94500 Champigny sur Marne, France - to justify the warranty.

### **I.3.1. Clause**

The following provisions are not exclusive to the benefit of the legal warranty provided by articles 1641 and following of the Civil Code, relative to hidden defects.

### **I.3.2. Purpose**

The warranty consists includes the free replacement of parts recognized as defective by our technical services.

Damage due to false manipulation, misuse, negligence or overload of the device, as well as those required due to customer's electrical network, power surges or faulty equipment, etc... are not covered by this warranty.

In addition, the use of products (liquids) not approved by the manufacturer makes this warranty null and void.

### **I.3.3. Transport**

The device must be transported in its original packaging and stored in a secure, dry place.

The device travels at risk of the user: in case of damage during transport, the recipient must file a claim with the carrier before accepting delivery of the device.

### **I.3.4. Duration**

This warranty remains in effect for 1 years from the date of invoice.

## I.4. Technical characteristics

Voltage :	<b>230V mono</b>	<b>115V mono</b>
Power:	1000W	1000W
Frequency :	50 à 60 Hertz	50 à 60 Hertz
Nominal intensity:	4,5A	9A

Weight: 5,9Kg  
 Package weight: 8,3Kg  
 Compatible bottle: 1 Litre  
 Maximum room volume to be treated: 1000m<sup>3</sup>  
 Average flow of the liquid: 1000 ml/h  
 Environmental conditions:

In condition of use :

- Temperature between 0°C and 40°C
- Ambient humidity between 10% and 90%
- Atmospheric pressure of 800 to 1060 hPa

Transport and storage conditions :

- Temperature between 0°C and 70°C
- Ambient humidity between 10% and 95%
- Atmospheric pressure of 700 to 1060 hPa

## II. PRECAUTIONS FOR USE

The manufacturer cannot be held responsible for any material or bodily harm resulting from the abnormal or faulty use of the device, or the disregard for guidance below:



Use only with NOCOTECH products.

- **Modifying the machine is forbidden.**
- To avoid risk of electric shock, this device must only be connected to a supply main with protective earth.
- Make sure that electric power meets regular standards.
- Hand-held radio frequency (use of remote control) communication devices may affect other electrical equipment and devices.
- This device may only use by professionals. According to requirements of norm NF EN 60601-1, it may cause radio interference and disturb operation of devices located close to the equipment.



- The use of accessories, transducers and parts, other than those supplied by OXY'PHARM for internal components, may result in increased emissions or decreased performance and robustness of the device.
- The device should not be stacked on or under other equipment. In case this is not possible, check the performance of the system in the desired configuration prior to treatment.
- Only OXY'PHARM products and liquids can be used with the OXY-30000.
- Do not introduce any object or liquid in the nozzle of OXY-30000.
- Check presence of nozzle and nozzle holder (cf. n°2 §III.1) prior to use.
- **Remove the bottle from OXY-30000 for transportation.**
- **Any intervention must be carried out by staff trained on the system.**

## II.1. Safety instructions

The device diffuses a chemical substance and makes noticeable sound, therefore it is necessary to pay attention to following recommendations.

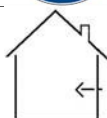
### **IMPORTANT RECOMMENDATIONS**



- Use only with NOCOTECH products.



- Always consult the accompanying documents.



- For indoor use only.



- Do not inhale the product.



- Do not direct the spray toward an incandescent source.



- Do not smoke close to the device.




- Grounded connections.



- No one should be in the room while the machine is operating
- Do not use in the presence of a flammable anesthetic mixture
- Handle the machine and bottle with care
- If a person must enter the room while the machine is running or during the required contact time, they must be equipped with individual protection equipment :
  - Noise-cancelling headphones.
  - Independent respirator.
  - Hydrogen peroxide resistant personal protective equipment.
- The individual should exit the room and seal it again as soon as possible.


# III. DESCRIPTION
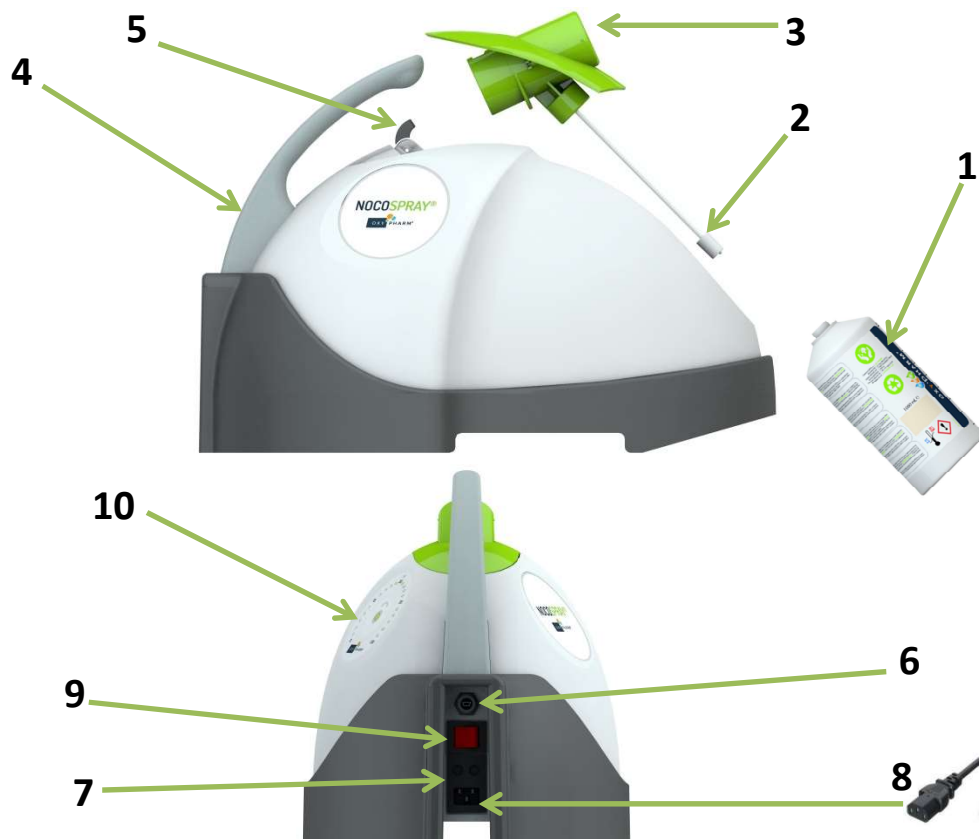
## III.1. OXY-30000 description



1	Bottle	6	USB port
2	Nozzle and nozzle holder	7	Fuses (5x20 10A T)
3	Nose of the diffuser (venturi)	8	Power cable
4	Handle	9	Power switch
5	Nozzle release button	10	Control touch pad

## III.2. OXY-30000 dimensions

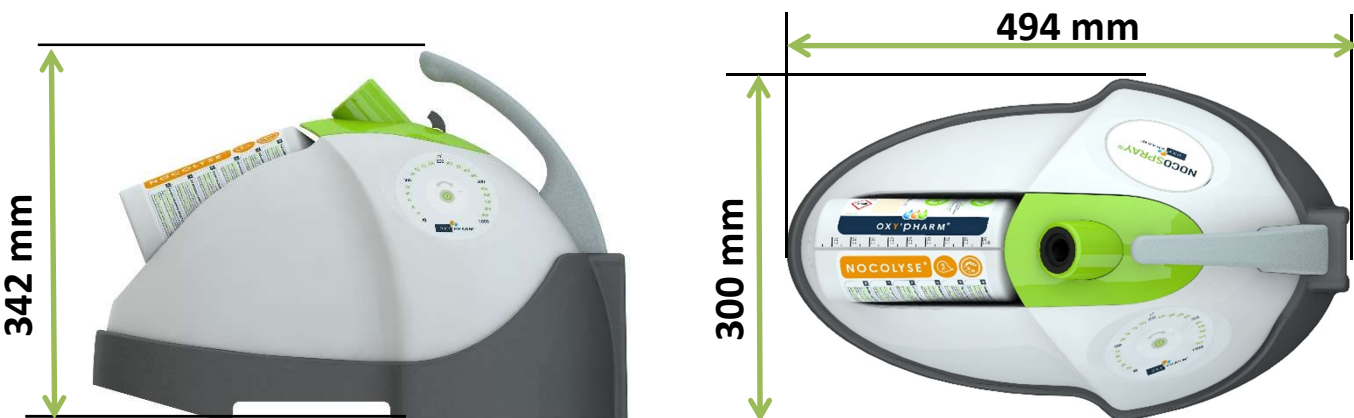

## IV. INSTRUCTIONS FOR USE

### IV.1. Positioning of the device

The OXY-30000 must be placed on a hard and flat surface and if possible in a corner of the room to treat.

#### IMPORTANT RECOMMENDATION






Before starting a cycle, check that there are no obstacles within 2 meters from the device.

### IV.2. Installation of the bottle



Follow instructions below to install the bottle.

		
Remove the diffuser nose by pressing the button	Open the bottle and screw on the nozzle (When opening a new bottle make sure to write the date in the area provided), do not over-tighten to avoid damaging the seal.	Place the nozzle back on the device. And press gently until you hear the « click ».

**Always check if the bottle contains enough product for the treatment  
(the device will not indicate whether there is sufficient liquid)**

## IV.3. Diffusion duration and consumption

Volume (m <sup>3</sup> )	Consumption (ml)	Duration of diffusion
10	10	36 seconds
15	15	54 seconds
20	20	1 minute & 12 seconds
25	25	1 minute & 30 seconds
30	30	1 minute & 48 seconds
35	35	2 minutes & 6 seconds
40	40	2 minutes & 24 seconds
45	45	2 minutes & 42 seconds
50	50	3 minutes
62.5	62.5	3 minutes & 45 seconds
75	75	4 minutes & 30 seconds
87.5	87.5	5 minutes & 15 seconds
100	100	6 minutes
112.5	112.5	6 minutes & 45 seconds
125	125	7 minutes & 30 seconds
137.5	137.5	8 minutes & 15 seconds
150	150	9 minutes
162.5	162.5	9 minutes & 45 seconds
175	175	10 minutes & 30 seconds
187.5	187.5	11 minutes & 15 seconds
200	200	12 minutes
225	225	13 minutes & 30 seconds
250	250	15 minutes
275	275	16 minutes & 30 seconds
300	300	18 minutes
325	325	19 minutes & 30 seconds
350	350	21 minutes
375	375	22 minutes & 30 seconds
400	400	24 minutes
425	425	25 minutes & 30 seconds
450	450	27 minutes
475	475	28 minutes & 30 seconds
500	500	30 minutes
525	525	31 minutes & 30 seconds
550	550	33 minutes
575	575	34 minutes & 30 seconds
600	600	36 minutes
650	650	39 minutes
700	700	42 minutes
750	750	45 minutes
800	800	48 minutes
850	850	51 minutes
900	900	54 minutes
950	950	57 minutes
1000	1000	1 hour



## IV.4. Functioning in standard mode

### IV.4.1. Use without remote control

EN



Connect the power cable and put the main switch on « I »  
(The light switches on)



Select the volume by dragging your finger on the clear part of the touchpad. (The range is between 10m<sup>3</sup> and 1000m<sup>3</sup>)



The light indicates the volume chose, if two lights are illuminated, the selected volume is the middle point between the two values.



After selecting the volume, press and hold the start button in the middle of the touchpad to start the diffusion cycle. A beeping sound will begin.



A beeping sound will continue for 15seconds and then the diffusion will begin. Leave and close off the room before the end of these 15 seconds.



The device can be stopped by pressing the start button again (device beeps when treatment ends).





## IV.4.2. Use with remote control

The remote control is used to start the device from outside the room, replacing the start button on control touch pad

### IV.4.2.1. Pairing with OXY-30000 (« blank » remote control)

When receiving the remote control and before first use, it needs to be paired with OXY-30000 appliance.

Pairing method :

		
<p>OXY-30000 needs to be turned off but plugged in</p>	<p>Press remote control button (+/- 5 seconds) remaining very close to the OXY-30000 device that needs to be paired until the Led flashes</p>	<p>While maintaining pressure on the remote control button, switch the OXY-30000 that needs to be paired on (both actions at the same time)</p>
		
<p>Remote control Led starts flashing quickly and the paired OXY-30000 beeps several times to signal that pairing was successful</p>		

#### IV.4.2.2. Functioning with the remote control

EN



Connect the power cable and put the main switch on « I »  
(The light switches on)



Select the volume by dragging your finger on the clear part of the touchpad. (The range is between  $10\text{m}^3$  and  $1000\text{m}^3$ )



The light indicates the volume chose, if two lights are illuminated, the selected volume is the middle point between the two values.



Press the remote control button to start the treatment.



A beeping sound will continue for 15 seconds and then the diffusion will begin.  
Leave and close off the room before the end of these 15 seconds.



The device can be stopped by pressing the start button again (device beeps when treatment ends).







#### IV.4.2.3. Pairing with another OXY-30000 (« re-pairing »)



The remote control can be paired to only one OXY-30000 at a time.

However, a remote control already paired can be 're-paired' to another OXY-30000

Re-pairing method :

		
<p>Check that all OXY-30000 close to the remote control are well switched off including the one to be re-paired</p>	<p>Plug in the OXY-30000 to be re-paired.</p>	<p>Press long on remote control button (+/- 30 secondes) remaining very close to the OXY-30000 device that needs to be re-paired until the Led flashes</p>
		
<p>During these 30 seconds, the Led light flashes and stops alternately. This phase is repeated 3 times</p>	<p>Once the Led light on the remote control flashes regularly and without stopping, switch the OXY-30000 that needs to be re-paired on, while maintaining pressure on the remote control button (both actions at the same time)</p>	<p>The Led light on the remote control starts flashing quickly and OXY-30000 beeps several times to signal that pairing was successful</p>

## IV.5. Functioning in timer mode

### IV.5.1. General informations

Plugged into OXY-30000 device and on mains socket, the timer enables postponing or repetition of treatments.  
Timer mode should be **selected to use OXY-30000 with a timer**.



Connect the power cable and put the main switch on « I »

To initiate timer mode, drag along the clear part of touchpad with your finger to the 1000m³ setting. Then continue in clockwise direction for 3 full rotations until three (3) lights are illuminated

Timer mode is on when 3 Led lights light up as follows :  
- Led 100m³ and 500m³ light up  
- Led 250m³ flashes.



Select the volume. As the selection is progressive, if two LEDs are lit, it is the intermediate value which is selected.  
**Note: selected volume appears only when the keyboard is touched/pressed**

Put the main switch on « 0 »

Plug the OXY-30000 socket into the timer and the timer in the mains socket



Set the timer, referring to the specific instructions for use.  
Nb : it can be set before it is plugged in

The timer will start and stop the OXY-30000 according to the setting

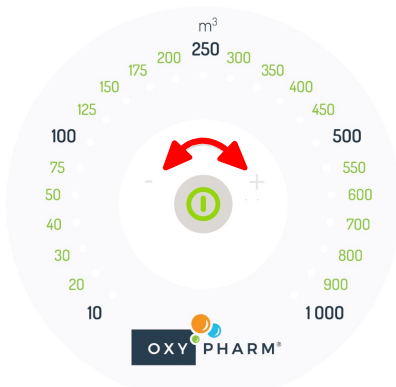
## IV.5.2. Timer mode features



If OXY-30000 is on timer mode when you switch it on, you can modify the volume or set the device back on standard mode **during 1 minute**.

**After this minute, OXY-30000 starts blinking for 15 seconds before starting treatment.**

When treatment is over, it is strongly advised to set the device back standard before switching it off.




When treatment is over, volume can be changed as long as the appliance is still switched on



The device cannot be stopped by pressing the start button on touchpad.

EN

## IV.5.3. Back on standard mode

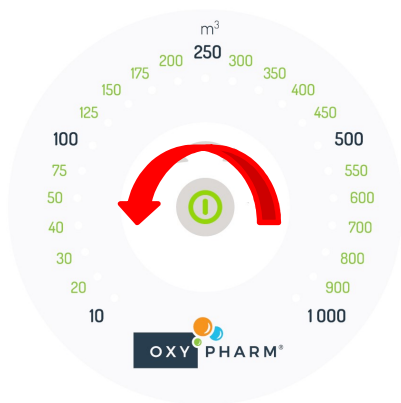

To set the appliance back on standard mode, touch the pad and drag your finger counterclock-wise back to the 10m3 position and then continue for an additional 3 full rotations

## V. TRACEABILITY SYSTEM

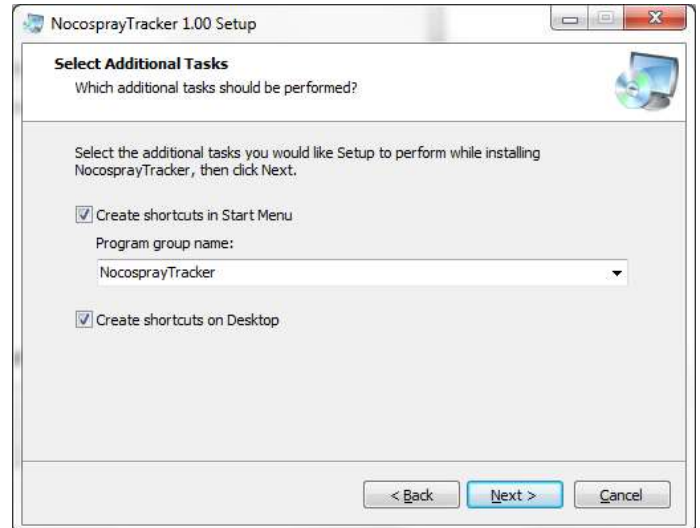
### V.1. Software installation

Click on “NocoSprayTacker\_Setup.exe” on the CD ROM or USB key supplied. Allow the program to install the software if Windows asks it.


1°) Click on « Next »



2°) Click again on « Next »




3°) Click on « Install »

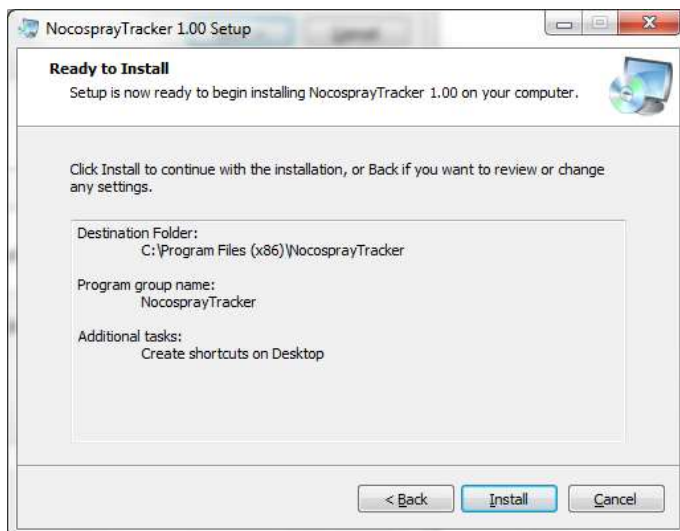

4°) Click on « Finish ».




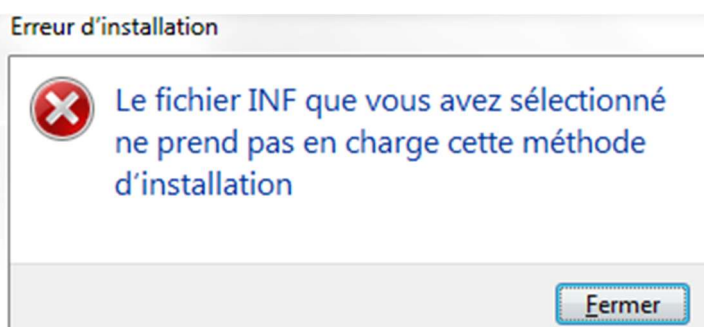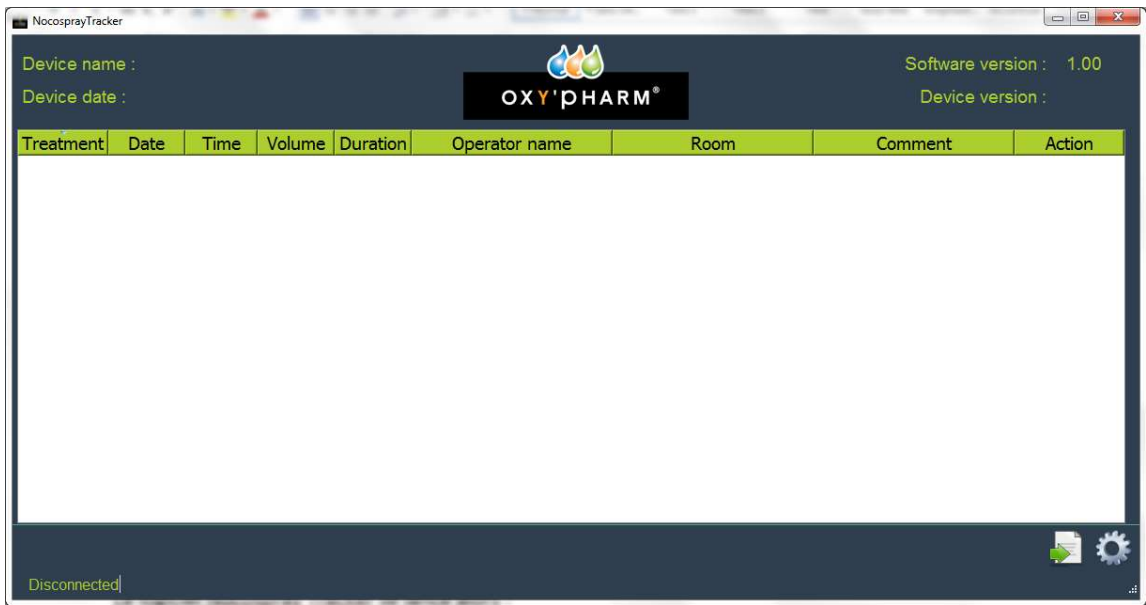If the following window appears, click on « Fermer »






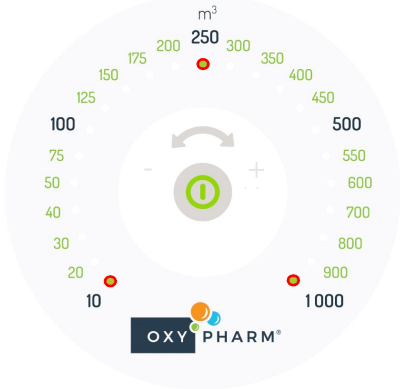

Nocospray Tracker software runs:



EN

## V.2. Transferring treatments historical

A mini USB connector enables transferring the treatments historical to a computer.

		
Plug the USB cable stick	The 10m3 and 1000m3 LEDs flash and the 250 m3 lights up.	<b>No treatment available in USB mode.</b> <b>To exit USB mode, unplug the cable from the computer</b>

Connect OXY-30000 on the USB port of your computer:

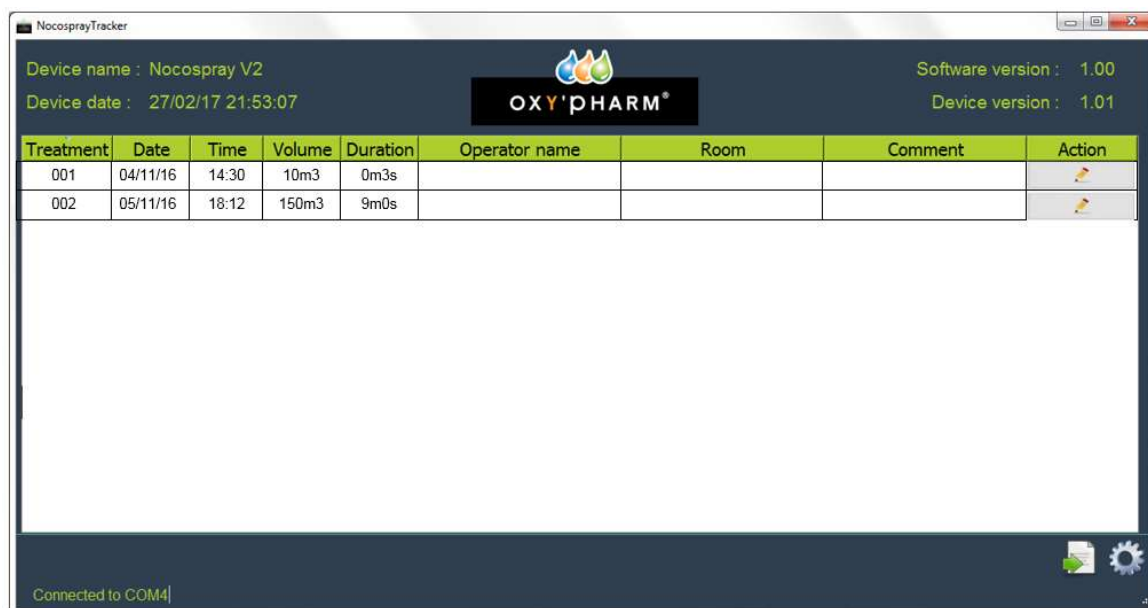
An icon on bottom right might appear to indicate that the device driver is running.



If your computer doesn't allow an automatic installation, the driver « .inf » is provided with the software, please contact your system manager.



Mention « Connected to COM4» on bottom left of the page, shows that Nocospray is well connected.



## V.3. READING AND MODIFYING DATA

### V.3.1.General information

Nocospray tracker software will launch automatically once installed. If it does not it can be initiated manually.

« Device Name »: Name of the device (editable parameter)

« Device date »: Date and time of the device (to be set the first time – it will automatically update afterwards)

« Software version »: Version of Nocospray Tracker software

« Device version »: Version of the intern software of Nocospray



Each line indicates each treatment. It is possible to organize each treatment following its number, date, time, volume, duration, room name, operator name or comment.

« Treatment »: treatment number

« Date »: Date of the treatment

« Time »: Hour of the treatment

« Volume »: Volume setting of the Nocospray

« Duration »: Duration of the treatment (this parameter is automatically calculated)

« Operator name »: Name of the operator who made the treatment.

« Room »: Name of the room where treatment has been launched

« Comment »: Other comment

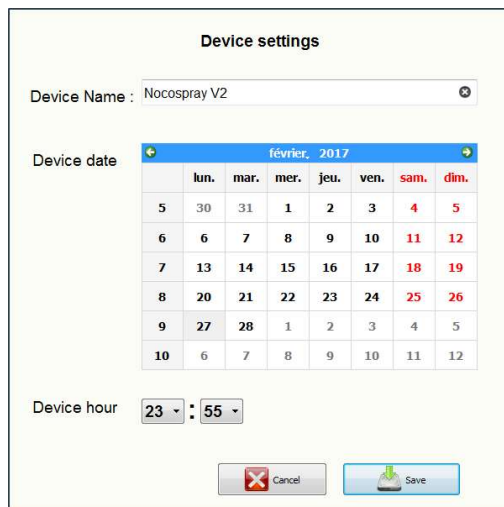


: Enter parameters « Operator name », « Room » and « Comment »

**WARNING: Once those parameters are saved, they cannot be modified**

### V.3.2.General settings modification

This icon gives you access to the following window:



**Device settings**

Device Name : Nocospray V2

Device date

	lun.	mar.	mer.	jeu.	ven.	sam.	dim.
5	30	31	1	2	3	4	5
6	6	7	8	9	10	11	12
7	13	14	15	16	17	18	19
8	20	21	22	23	24	25	26
9	27	28	1	2	3	4	5
10	6	7	8	9	10	11	12

Device hour 23 : 55

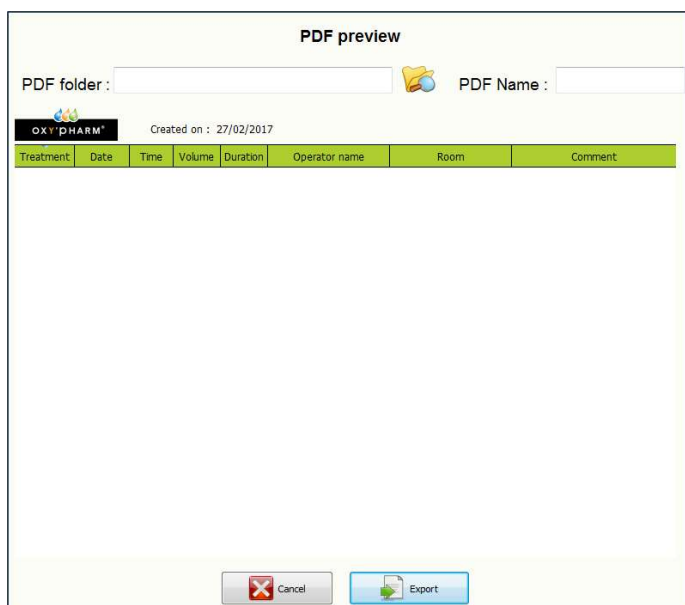
Cancel Save

On this window device name « Device Name », date and hour of the device can be modified.


### V.3.3.Export records (treatment historical)



This icon enables exporting the all recorded treatments to a PDF file as long as they have not been deleted. It brings you to the following window :




**PDF preview**

PDF folder :   PDF Name :

OXI-PHARM® Created on : 27/02/2017

Treatment	Date	Time	Volume	Duration	Operator name	Room	Comment
-----------	------	------	--------	----------	---------------	------	---------

Cancel Export

Using this icon  , you can choose the folder where you want save the data (PDF file).

« PDF Name » : enter the name of the record file.

Then, click to « export » to record file.

### V.3.4. USB mode exit

To exit USB mode, unplug the cable from the computer.

## VI. TROUBLESHOOTING

The table below shows a list of common problems and likely solutions. If the problem persists, contact your reseller or our Customer Service Center.



**Before any intervention, switch off the device and disconnect the power cord**

Failure	Resolution
The device does not switch on	<ul style="list-style-type: none"><li>- Check the power cable connection</li><li>- Check if the switch is in « I » position and if the light is on</li><li>- Check fuses</li></ul>
The device works but the consumption is incorrect	<ul style="list-style-type: none"><li>- Check that nozzle is attached securely, is in working, and not blocked</li><li>- Check if the bottle contains enough liquid</li></ul>
The device has fell and has a malfunction	Contact the Customer Service Center

## VII. HYGIENE AND MAINTENANCE


The OXY-30000 device requires regular maintenance.



- **Switch off the device before proceeding any cleaning.**
- **Do not use abrasive products, bleach, acetone or solvents for maintenance or cleaning operations of the device. These aggressive agents may damage the surface materials (discoloration, tasks, drips, microcracks...)**
- **For the cleaning of external surfaces, use a slightly damp cloth or wipes.**

## VIII. EMC COMPLIANCE BOARD

EMC Compliance according to IEC 60601-1-2 (2007)			
The <b>OXY-30000</b> is intended for use in the electromagnetic environment specified below. The customer or the user of the <b>OXY-30000</b> should assure that it is used in such an environment.			
Emission test	Standard	Compliance	Electromagnetic environment
Emissions RF	CISPR 11	Group 1	The <b>OXY-30000</b> uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.  The <b>OXY-30000</b> is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Emissions RF	CISPR 11	Class B	
Harmonics emissions	IEC 61000-3-2	N/A	
Voltage/Flicker test	IEC 61000-3-3	N/A	
The <b>OXY-30000</b> is intended for use in the electromagnetic environment specified below. The customer or the user of the <b>OXY-30000</b> should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment
Electrostatic discharges IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical Fast Transients IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle 40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles 70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 s	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle 40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles 70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 s	Mains power quality should be that of a typical commercial or hospital environment. If the user of the system requires continued operation during power mains interruptions, it is recommended that the system be powered from an uninterruptible power supply or a battery. $U_T$ is the a.c. mains voltage prior to application of the test level.
Magnetic field at power frequency IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Conducted RF IEC 61000-4-6	3 Vrms 150kHz-80MHz	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the NOCOSPRAY 2, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance : $d = 1,67 \cdot \sqrt{P}$
Radiated RF IEC 61000-4-3	3V/m 80MHz-2.5GHz	3V/m	$d = 1,67 \cdot \sqrt{P} \text{ 80MHz-800MHz}$ $d = 2,33 \cdot \sqrt{P} \text{ 800MHz-2,5GHz}$ <p>where <math>P</math> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <math>d</math> is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey (a) should be less than the compliance level in each frequency range (b)</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
(a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the NOCOSPRAY 2 is used exceeds the applicable RF compliance level above, the NOCOSPRAY should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the NOCOSPRAY 2.			
(b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the OXY-30000.			
The <b>OXY-30000</b> is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the <b>OXY-30000</b> can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the <b>OXY-30000</b> as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150kHz to 80MHz	80MHz to 800MHz	800MHz to 2.5GHz
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.737
1	1.167	1.167	2.330
10	3.690	3.690	7.368
100	11.67	11.67	23.300
For transmitters rated at a maximum output power not listed above, the recommended separation distance $d$ in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

## IX. RECYCLING

### IMPORTANT INFORMATION ABOUT THE ELIMINATION COMPATIBLE WITH THE ENVIRONMENT

This device is subject to separate collection for electrical and electronic equipment in accordance with Directive 2002/96 / EC (applicable in European Union countries and other European countries with separate collection systems).

This symbol that appears on the device indicates that at the end of its life cycle, the device will need to be treated separately from household waste.

The device will be assigned to a center of selective collection for electrical and electronic appliances, or reported to the dealer on the purchase of a new device.



The user is responsible for the delivery of the device at the end of life to the collection authorities or to a corporation duly empowered under penalty of sanctions provided for by the legislation on the elimination of waste. The waste supplies must be managed in compliance with the local regulations relating to the management of waste.

Proper differential collection, and the subsequent recycling, processing and environmentally compatible disposal of waste equipment avoids unnecessary damage to the environment and possible related health risks and also promotes recycling of the materials used in the appliance

## BELOW FOR CHINESE MARKET ONLY

WE STRONGLY recommend customer in China contact our legal distributor of OXYPHARM in China: Shenzhen RunLian Eco-Tech for more technical information and after-sales services.

**OXY-30000** is our latest dry-fogging system with most advanced fogging technology and more convenience.

**TO ENSURE YOUR BENEFITS, YOU MAY ASK LOCAL DISTRIBUTOR PROVIDE:**

**-LEGAL SN NUMBER**

**-LEGAL CERTIFICATE OF ORIGIN ( TYPE CLEARLY)**

我们强烈建议我们在中国的客户联系我们在中国的合法代理：深圳润联环保科技有限公司 以获取更多技术支持和售后服务

OXY-30000是我们最新的干雾系列产品，凝聚着最先进的雾化技术和更多的便利性

为确保您的利益，您可以向当地经销商索取：

- 合法的SN号码
- 合法的原产地证明

中国代理商售后及技术支持热线：

400—6787—918 （ 24 hours）

0755—29105906 （8：30—18：00 working day）

或扫描右方二维码：

